serenysens

Manuel d'utilisation
Users' manual

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TO REDUCE THE RISK OF ELECTRIC SHOCK, FIRE AND INJURY, THIS PRODUCT SHOULD BE USED IN ACCORDANCE WITH THE FOLLOWING INSTRUCTIONS :

This appliance may be used by children from the age of 8 and by persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge, provided that they are under supervision or have been instructed in the safe use of the appliance and understand the hazards involved. Children must not play with the appliance. Cleaning and maintenance by the user should not be performed by children without supervision.

CAUTION : Do not use this appliance near a bathtub, shower, washbowl or other receptacle containing water.

- This device is intended for home use only.
- This device is not a toy.
- Check the device before and after each use. Do not use it if the product does not work properly.
- Do not use this product if any parts are missing, damaged or worn.
- Stop the device immediately, if you feel discomfort, pain or if you notice something abnormal (smoke, odors or particular sounds).
- In order to prevent any allergies (redness, skin reaction) please test the belt on your skin for 30 minutes. If nothing appears, you can use your equipment freely. Otherwise, please remove the belt immediately and

consult a doctor. Failure to follow these contraindications could put you and your health at risk. If you have any other concerns about the heating belt and your health, contact your physician before your first use.

- Do not leave the device unattended when it is plugged in.
- Do not handle or store the device near flammable materials.
- Keep away from all sources of heat, including direct sunlight.
- Do not store the unit in an environment that is too cold or too humid to prevent damage.
- Do not store the unit in an environment exposed to direct sunlight.
- Store the unit in a dry, dust-free place.
- Keep the unit clean, dust and dirt free.
- Do not immerse the device in water. Do not handle with wet hands.
- The unit should not be used for heating animals.
- Make sure that the product and the power supply are completely dry before using the product.
- Do not use it improperly. Make sure it is not knocked or dropped.
- If the product is damaged, contact the store, its service department, or a person of similar qualification.
- If the power cord is damaged, it should be replaced by the manufacturer, its service agent or similarly qualified persons.
- Do not attempt to repair or modify the unit or any of its components yourself. Contact the store, its service agent or a person of similar qualification.
- The unit must be powered according to the recommendations in the manual and on the product label.
- Do not unplug the machine by pulling on the cord, pull on the plug.
- Always disconnect this unit from the mains immediately after use and before cleaning or maintenance.
- People with sensitive skin, pregnant women and people with scratches or cuts should not use it.
- It should not be used for a long period of time (maximum 90 minutes).
- It must be unplugged when not in use or before going to sleep.
- Never place the product on leather materials such as a sofa or metal materials when it is switched on or hot.
- Never place the remote control on or in the product when it is turned on or hot (risk of overheating).
- Children under the age of three are not to use this appliance due to their inability to react to overheating
- The appliance cannot be used by children older than 3 years and younger than 8 years due to their inability to react to overheating
- If the appliance is to be used continuously, the interval should be more than two hours
- The appliance is only to be used with the power cord FY-201 which is supplied in this pack.
- When storing the appliance, allow it to cool down before folding
- Do not crease the appliance by placing items on top of it during storage
- Examine the appliance frequently for signs of wear or damage. If there are such signs, if the appliance has been misused or does not work, return it to the supplier before switching it on again
- This appliance is not intended for medical use in hospitals
- Do not use if wet
- This appliance must not be used by persons insensitive to heat and other very vulnerable persons who are unable to react to overheating.

Contraindications of the lumbar heating belt :

In case of doubt and particular health problems, it is recommended to seek the advice of a health professional.

Contraindicated in all cases of :

- Pregnant women.
- Stimulation of the thoracic area if the person has a cardiac arrhythmia.
- In case of abdominal or inguinal hernia.
- Autoimmune diseases.
- Metabolic diseases.
- Cardiovascular diseases.
- Respiratory diseases.
- Taking medication or food supplements.
- In case of fever
- In case of open skin wounds and eczema
- Serious illnesses, bacterial infections or inflammatory processes
- Sun or general burns
- Skin intolerances
- Recent operations
- Thrombosis, thrombophlebitis
- Stents that have been active for less than six months
- Arteriosclerosis, arterial circulatory disorders
- Untreated arterial hypertension
- Cardiac arrhythmias
- Other cardiovascular diseases
- Bleeding disorders (hemophilia)
- Severe forms of diabetes mellitus
- Cancer or tumor diseases
- Acute arthritis
- Neurological diseases and epilepsy

- Progressive muscular dystrophy
- Abdominal or inguinal hernias
- Undiagnosed back pain
- Diseases of internal organs

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- This package contains:
 - 1 heating belt
 - 1 cable with integrated remote control and power plug
- Input : 220-240V~50/60Hz
- Power : 100W
- Lifetime : 2000H
- Number of switching cycles : 40 000 fois
- Cable length : 2.4m
- Material : 100% Polyester, ABS plastic and PVC
- Dimensions:
 - Adjustable belt from 65 to 130 cm, height : 28 cm
- Weight : 390grs

NETTOYAGE ET STOCKAGE / CLEANING AND STORAGE

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- Always unplug the power cord when the belt is not in use.
- Turn off the belt before cleaning.

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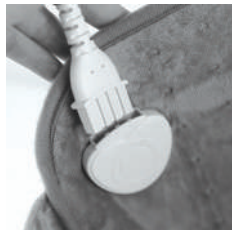
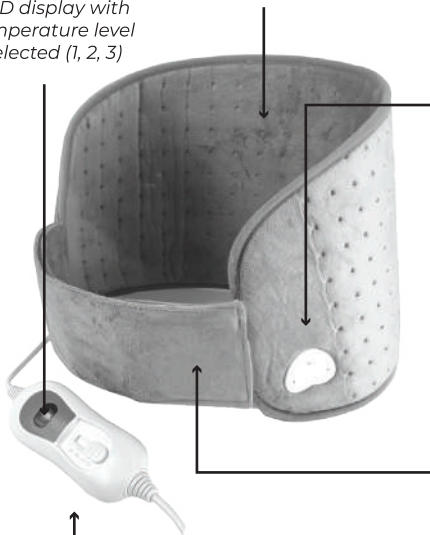
1. Connect the power cable to the belt and the other end to a power outlet.
2. Fasten the belt around your waist with the self-gripping straps so that the heating part is in the back.
3. Now sit back and choose your temperature (between 1: 50°, 2: 55° and 3: 60°), using the remote control.

4. When you are finished using it, turn the button to "0".

- Remember to turn the switch back to the "0" position if you are not using it anymore.
- If you want to turn it on again, turn the button to "0", wait 10 seconds and then turn the belt on again.

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EN: If you wish to dispose of the packaging, please observe the relevant environmental regulations in your country.

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EN : Electrical products should not be disposed of with household waste. Please recycle them at the designated collection points. Contact your local authority or dealer for advice on recycling.

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Keep the receipt. If your product does not work properly, has manufacturing or material defects, please contact the shop. The warranty applies to material or manufacturing defects. Wear parts that are subject to normal wear and tear are excluded from the warranty. This warranty is void if the device has not been used or maintained properly. For proper use of the product, all instructions in the operating instructions must be followed precisely. The product is intended exclusively for private, non-commercial use. Unauthorized and improper handling, use of force and interventions not carried out by our authorized service center will result in loss of warranty.

serenysens

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Copyright © Reproduction interdite. MADE IN P.R.C. - FABRIQUÉ EN R.P.C.